


## Instructions

# Motor, Heat/Massage


**Haley Recliner / Hannah Patient Chair**

### Electrical Connections

HPFi

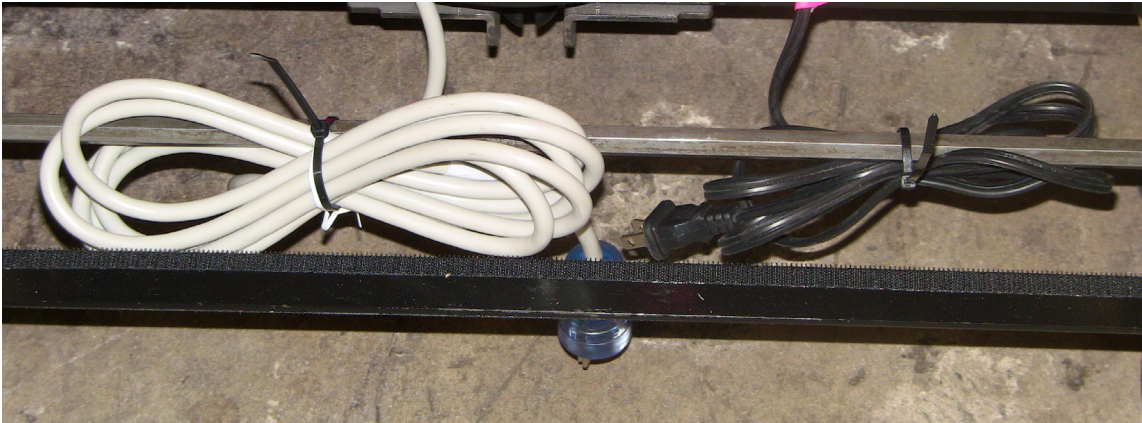

Step 1: Carefully cut the cable tie securing the AC power cable from the cross bar. Make certain that the power cords are firmly engaged. The cord for Heat & Massage labeled A should be plugged in to the power module also labeled A. The cord for the Haley Comet motor labeled B should be plugged in to the power module labeled B.

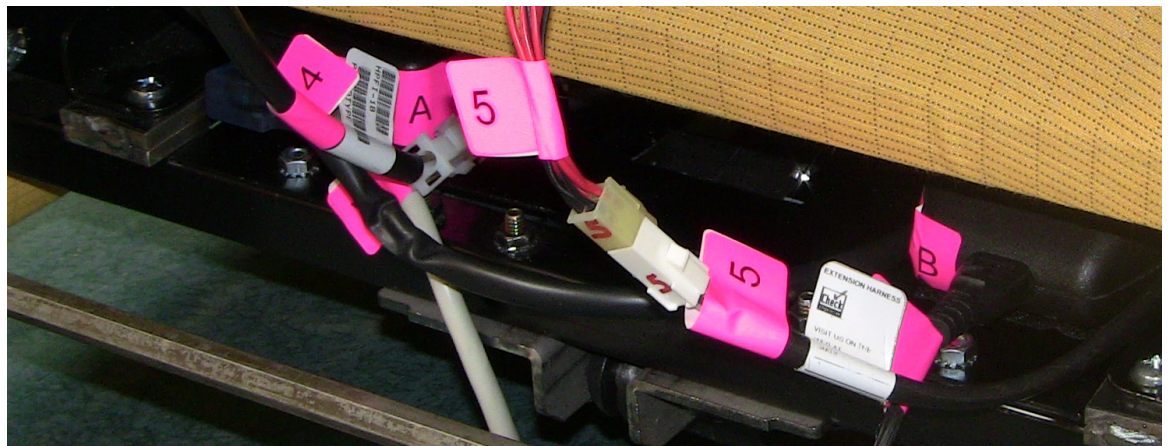

Step 2: Connect cable assembly harnesses from the seat to the back making certain that the numbers correspond (5 connects to 5, 4 connects to 4.)